

Medical  
Supporting  
Equipment

# Medicop Medical Equipment

Dear clients,

With present brochure we would like to introduce our newest medical rail system, which comply with demands of relevant international standards and include the most advanced technical solutions, based on our 30 years experiences with production of medical equipment. In this brochure you will find the most important benefits which the system offers, while detailed technical data and drawings are available in our Technical Catalogue.

Our entire production program of medical equipment includes the following product-groups:

Central System  
of Medical Gases

Medical  
Suction Units

Oxygen  
Therapy Devices

Bed Head  
Units

Medical  
Supporting Equipment

## Content


4 MEDICAL  
RAIL

6 SHELVES FOR  
RAIL SYSTEM

8 INFUSION  
STANDS

10 PARTITION  
CURTAINS

12 BASKETS

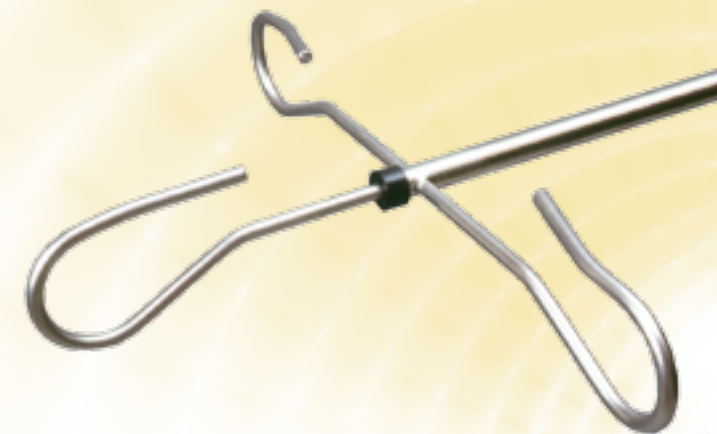

# MEDICAL RAIL

Made of stainless steel

## MEDICAL RAIL


Rail system is used for hanging or attaching different medical instruments to the wall rail. The rail of 25x10 mm cross section can be attached directly to the wall or on a pendant. The rail is made of stainless steel and includes connections for potential equalization.

The length of the rail can be specified by the customer. On the both ends the rail is closed by appropriate covers.

TECHNICAL DATA	
<b>Material:</b>	stainless steel
<b>Cross section:</b>	25 x 10 mm
<b>Load capacity Qmax:</b>	100 kg/m
<b>Recommended mounting height:</b>	lower rail: 50 cm upper rail: 130 cm

AVAILABLE MODELS	
<b>1600003</b>	Wall rail made of stainless profile, 25 x 10 mm, including connecting accessories
<b>1600001</b>	Wall rail made of aluminium profile, 25 x 10 mm, including connecting accessories

**1600003** Wall rail made of stainless profile, 25 x 10 mm, including connecting accessories all tubes, filters and connectors

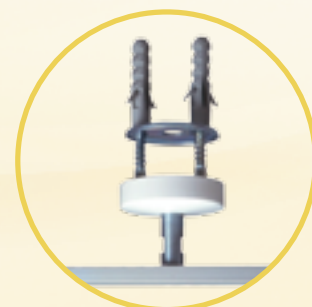

## SHELVES FOR RAIL SYSTEM

The shelves are intended for depositing EKG monitors, defibrillators and other medical devices, as well as various documents.

The shelves are fixed to the rail with two accessories and could be additionally supported beneath with screws. Standard dimensions of the shelves for monitors are 530x370 mm or 370x420 mm, while standard dimension for the shelves for documents is 330x260 mm.

Customer specified shelf dimensions are also available. All shelves are made of stainless sheet-metal.

TECHNICAL DATA	
<b>Material - shelves:</b>	stainless steel
<b>Material - drawer:</b>	ABS
<b>Load capacity:</b>	20 kg (models 1610001, 1610002)

AVAILABLE MODELS	
<b>1610001</b>	Shelf for monitor with drawer; 530 x 370 mm
<b>1610002</b>	Shelf for monitor; 530 x 370 mm
<b>1610005</b>	Shelf for documents; 330 x 260 mm
<b>1610027</b>	Vertical rail carrier of documents 220 x 300 mm



**1610001** Shelf for monitor with drawer; 530 x 370 mm



**1610002** Shelf for monitor; 530 x 370 mm



**1610005** Shelf for documents; 330 x 260 mm




**1610027** Vertical rail carrier of documents 220 x 300 mm

## BASKETS


Baskets are made of chromium plated net like sheet metal and are intended for depositing different sanitary material, especially all types of catheters.

They can be hung on the rail or any other carrier. Standard basket dimensions are: 80x140x470 mm, 200x250x300 mm and 100x200x250 mm. Other dimensions of baskets are available upon request.


AVAILABLE MODELS	
<b>1610040</b>	Basket for catheters for wall rail (80 x 140 x 470 mm)
<b>1610041</b>	Basket for catheters for wall rail (200 x 250 x 300 mm)
<b>1610042</b>	Basket for catheters for wall rail (100 x 200 x 250 mm)
<b>Remark:</b>	Other dimensions of baskets are available upon request.

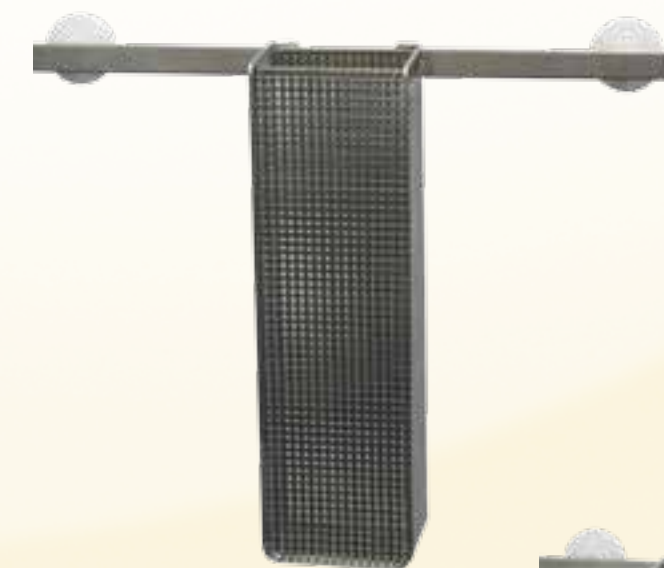

**1610040** Basket for catheters for wall rail (80 x 140 x 470 mm)

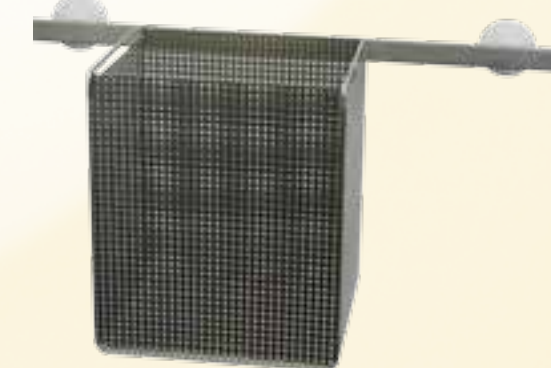

**1610041** Basket for catheters for wall rail (200 x 250 x 300 mm)

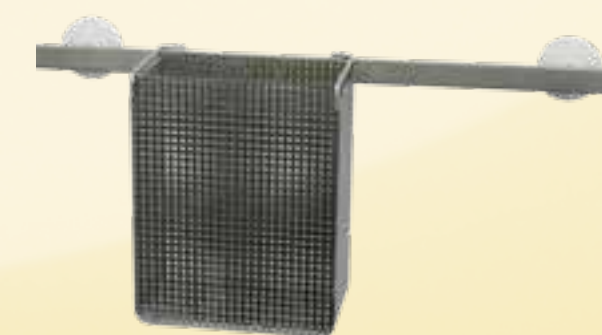

**1610042** Basket for catheters for wall rail (100 x 200 x 250 mm)

# MEDICAL RAIL

Made of stainless steel

## PARTITION CURTAINS

The carrier of curtain can be fixed to the wall-rail or to the wall directly. The pull out telescopic carrier is in folded position 700 mm long, while the pulled-out length is 1950 mm.

The height of the partition curtain is 1400 mm and its standard color is white. Curtains of other colors are available upon request.

TECHNICAL DATA	
<b>Material of telescopic carriers:</b>	Eloxized aluminum bars
<b>Dimensions – folded position:</b>	700 mm
<b>Dimensions – pulled out length:</b>	up to 1950 mm
<b>Height of the partition curtain:</b>	1400 mm
<b>Standard colour of curtain:</b>	white

AVAILABLE MODELS	
<b>1610035</b>	Carrier of a partition curtain with a curtain, rail mounted
<b>1610036</b>	Carrier of a partition curtain with a curtain, wall mounted
<b>1610037</b>	Carrier of a partition curtain without a curtain, rail mounted
<b>1610038</b>	Carrier of a partition curtain without a curtain, wall mounted
<b>1610015</b>	Carrier of disposable bags made of stainless steel, with wall rail bracket
<b>1610030</b>	Folding rib tube holder, for wall rail
<b>Remark:</b>	Curtains of other colours are available upon request.




## INFUSION STANDS - RAIL-MOUNTED OR MOBILE

Infusion stands serve for hanging infusion liquid jars and infusion pumps. Stands can be rail-mounted or mobile - with a trolley. Rail mounted stands have three hooks for hanging of infusion bottles, whereas the mobile one has four hooks.

The height of rail mounted infusion stands is up to 1200 mm above the rail, and the height of the mobile stand is 1250 to 1950 mm.


The height is adjustable by screw or by button. Infusion stands are made of stainless steel bars.

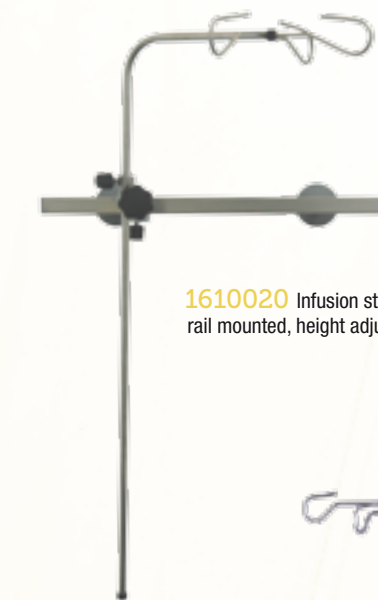

**1610020** Infusion stand with 3 hooks, rail mounted, height adjustable by screw

TECHNICAL DATA	
<b>Material:</b>	stainless steel
<b>Height - rail mounted infusion stands</b>	up to 1200mm above the rail
<b>Height of the mobile stand:</b>	1250 to 1950 mm

AVAILABLE MODELS	
<b>1610020</b>	Infusion stand with 3 hooks, rail mounted, height adjustable by screw
<b>1610021</b>	Infusion stand with 3 hooks, rail mounted, height adjustable by button
<b>1610022</b>	Mobile infusion stand with 4 hooks, including trolley with 5 antistatic wheels

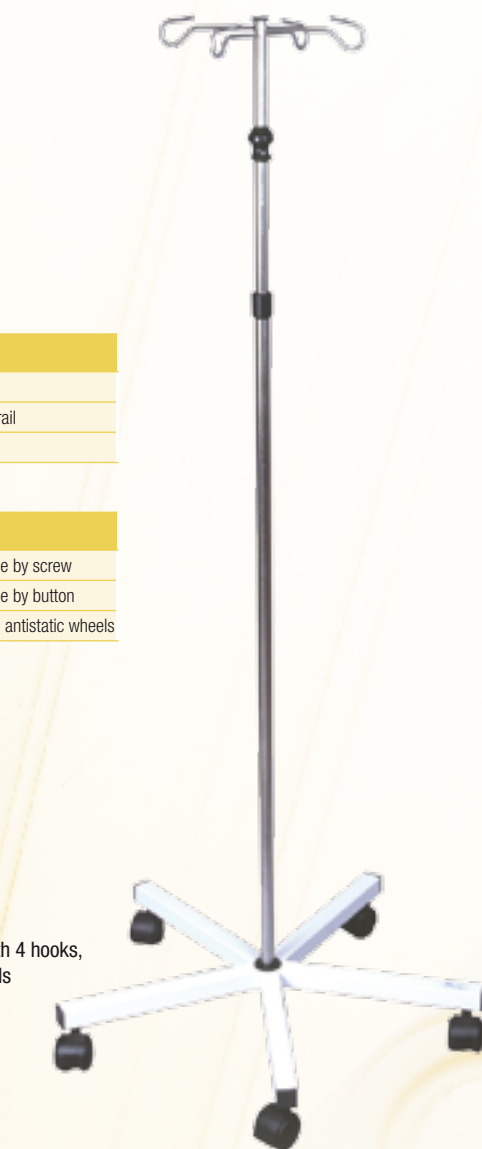

**1610022** Mobile infusion stand with 4 hooks, including trolley with 5 antistatic wheels



MEDICOP medicinska oprema d.o.o.  
Obrtna 43 (p.p. 161)  
SI - 9000 Murska Sobota, Slovenia  
T: +386 2 53 91 250 | F: +386 2 53 91 255  
info@medicop.si | www.medicop.si